


The background of the slide is a light gray gradient with several realistic water droplets of various sizes scattered across it. The droplets have highlights and shadows, giving them a three-dimensional appearance. The text is centered in a bold, red, sans-serif font.

FOREARM X RAY
BY
S.S.SIVA SAGAR
TEACHING ASSISTANT IN
RADIOLOGY
CENTURIAN UNIVERSITY

VIEWS

- FOREARM AP.
- FOREARM LATERAL.

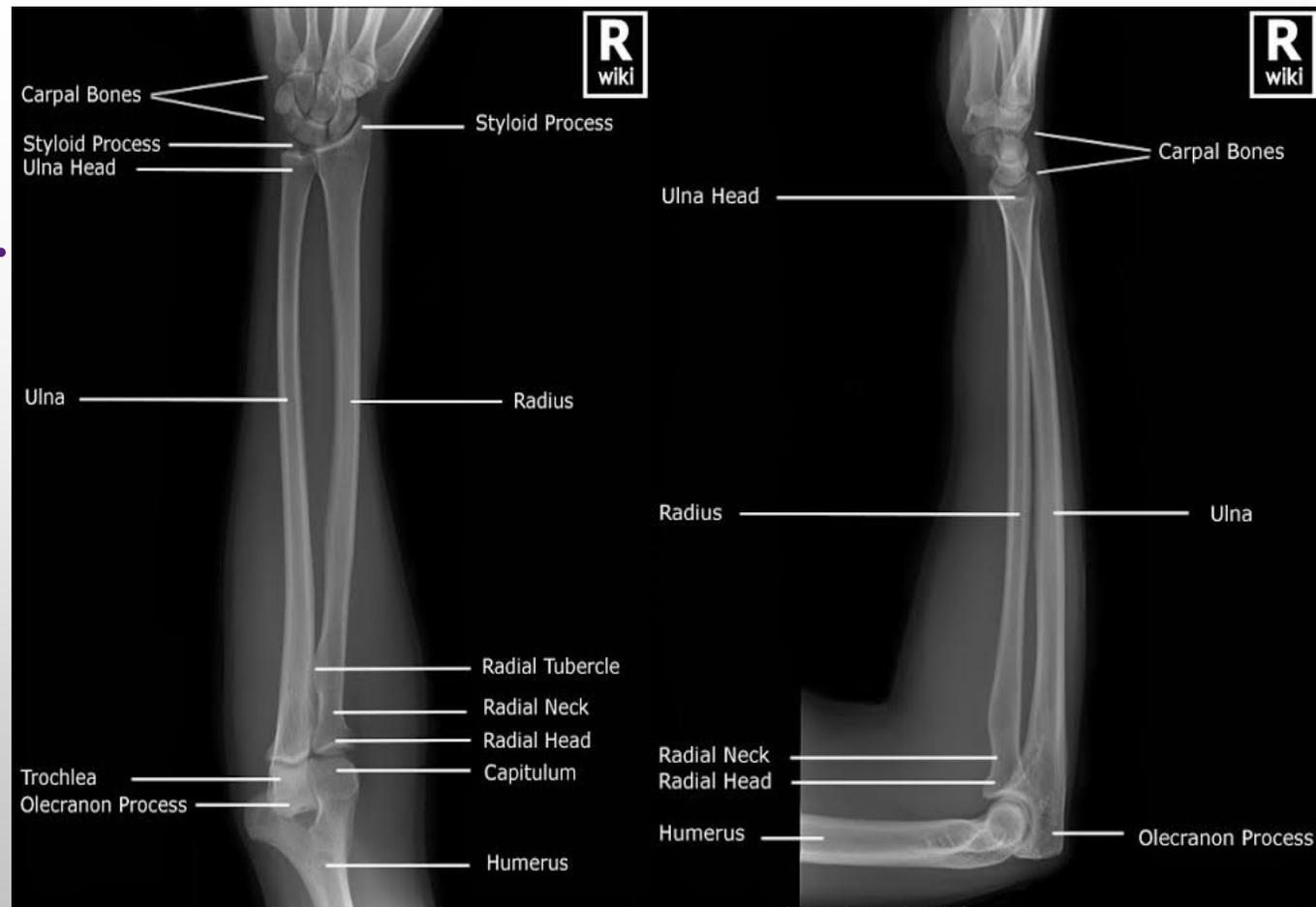

INDICATIONS

- **PAIN.**
- **TENDERNESS.**
- **SWELLING AND DEFORMITY OF FOREARM.**
- **SUSPECTED FRACTURES.**
- **POST OP CASES.**
- **SUSPECTED FOREIGN BODY.**



CONTRAINDICATIONS

- **PREGNANCY**
- 

PREPARATION

- **COLLECT HISTORY FROM PATIENT.**
- **PATIENT IS ASKED TO REMOVE ALL RADIOPAQUE MATERIALS.**
- **WEAR HOSPITAL GOWN.**
- **IF FEMALE PATIENT COME TO X RAY , CHAPERONE IS NECESSARAY.**
- **SHIELDING IS PROVIDED TO PATIENT.**


FOREARM AP POSITIONING

- **THE PATIENT IS SEATED ON THE SIDE OF X RAY TABLE.**
- **THE FOREARM IS FULLY SUPINATED, ELBOW EXTENDED AND PLACED ON THE CASSETTE.**
- **THE HAND IS IMMOBILISED WITH SAND BAG, SIDE IS MARKED AND BEAM IS COLLIMATED.**



FOREARM LATERAL POSITIONING

- **PATIENT IS SEATED ON THE SIDE OF X RAY TABLE.**
- **THE ELBOW FLEXED 90 DEGREES AND IT IS PLACED LATERALLY ON CASSETTE.**
- **THE HAND IS IMMOBALISED BY SAND BAGS, SIDE IS MARKED AND BEAM IS COLLIMTED.**

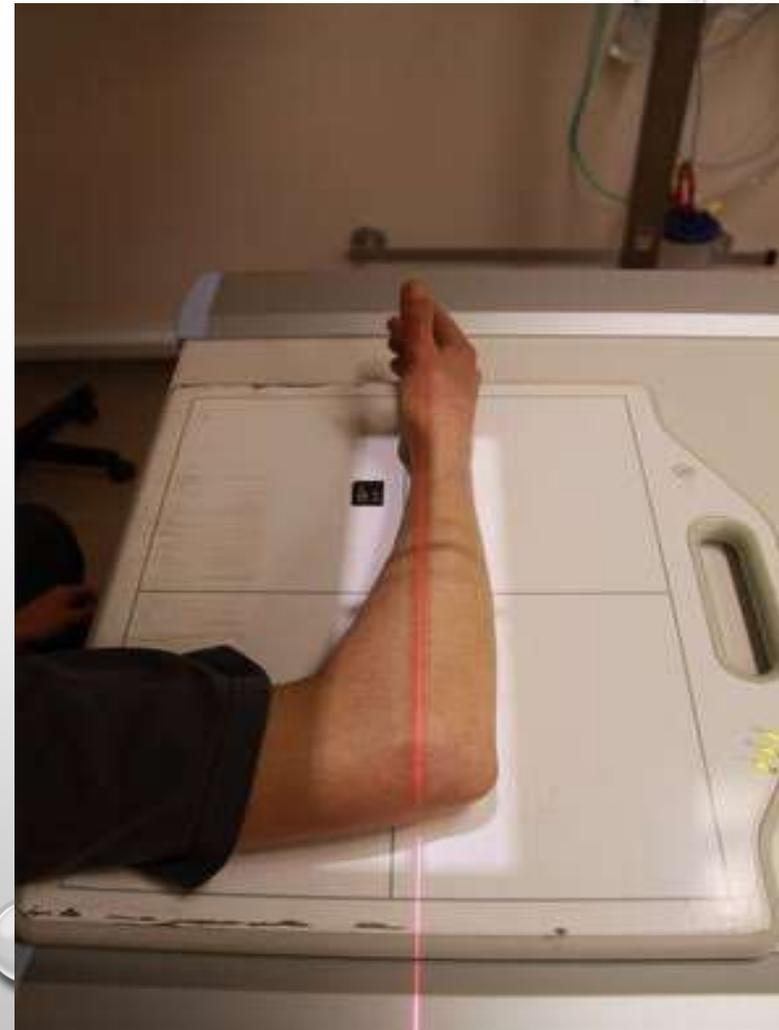


CENTERING

- ***FOREARM AP : IT IS DONE ON ANTERIOR ASPECT OF MIDDLE OF FOREARM.***
- ***FOREARM LATERAL : IT IS DONE ON THE LATERAL ASPECT OF MIDDLE OF FOREARM.***

TECHNIQUE

- **KVP : 80**
- **MAS : 50**
- **FSS SMALL**
- **FFDD : 100 CMS**
- **WITH OUT GRID.**
- **CASETTE SIZE : 12" X 10"**

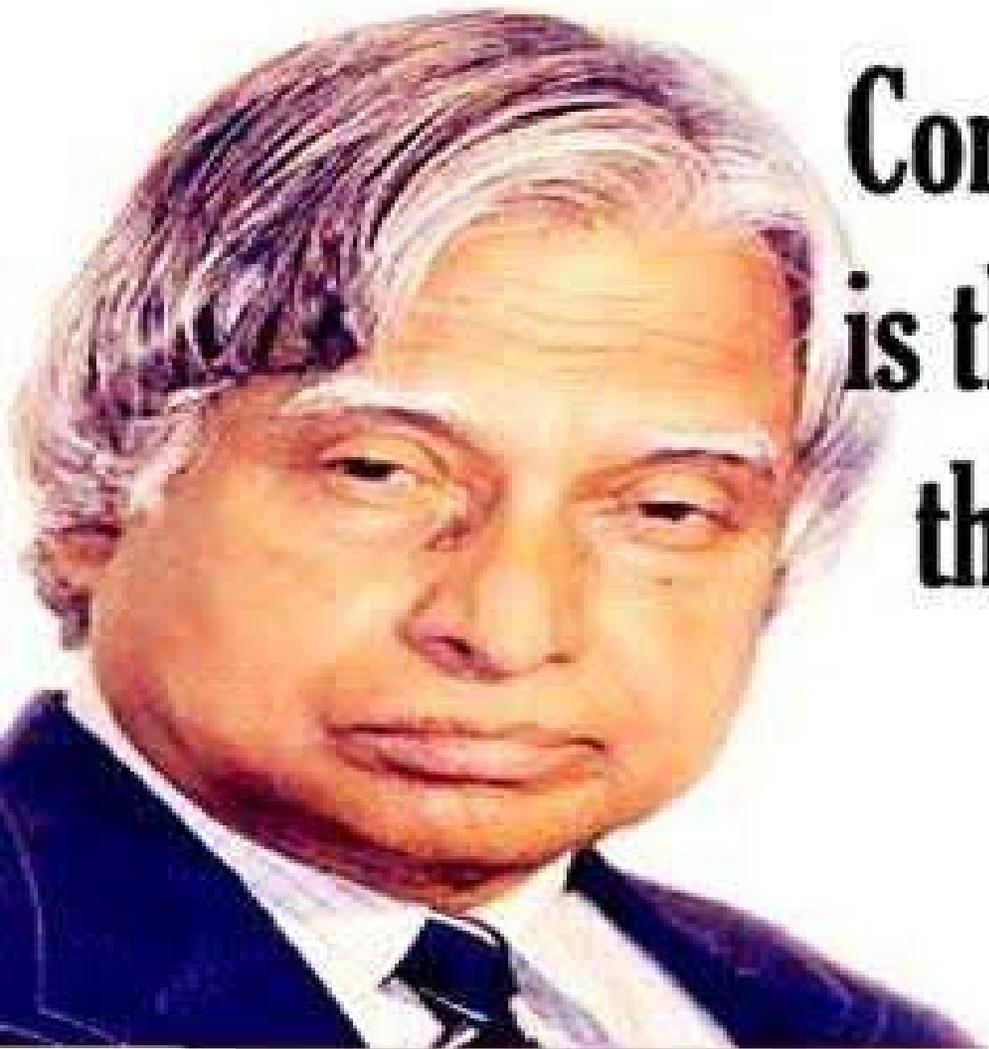

**Confidence and Hard-work
is the best medicine to kill
the disease called failure.**

**It will make u a
successful person**